

COBAS AMPLISCREEN HBV TEST

COBAS AmpliScreen™ HBV Test	HBV	96 Tests	P/N: 21118315018
COBAS AmpliScreen Multiprep Specimen Preparation and Control Kit	MULTIPREP/ CTL	96 Tests	P/N: 3302555018
COBAS AMPLICOR™ Wash Buffer	WB	500 Tests	P/N: 20759899123

TABLE OF CONTENTS

INTENDED USE.....	3
SUMMARY AND EXPLANATION OF THE TEST.....	4
PRINCIPLES OF THE PROCEDURE	5
Specimen Processing	5
PCR Amplification.....	6
Selective Amplification	7
Hybridization Reaction	7
Detection Reaction.....	8
MATERIALS PROVIDED BY ROCHE	8
COBAS AmpliScreen Multiprep Specimen Preparation and Control Kit.....	8
COBAS AmpliScreen HBV Test.....	9
COBAS AMPLICOR Wash Buffer Kit.....	9
OTHER MATERIALS REQUIRED BUT SOLD SEPARATELY (MAY BE PURCHASED FROM ROCHE)	10
MATERIALS REQUIRED BUT NOT PROVIDED BY ROCHE	11
REAGENTS.....	12
STORAGE INSTRUCTIONS	16
PRECAUTIONS	17
REAGENT PREPARATION	19
SPECIMEN COLLECTION, STORAGE AND POOLING	22

PROCEDURAL NOTES	24
INSTRUCTIONS FOR USE	28
Amplification, Hybridization and Detection of Stored Processed Specimens	28
Hybridization and Detection of Stored Denatured Amplicon.....	28
QUALITY CONTROL PROCEDURES	37
INTERPRETATION OF RESULTS	38
Invalid Test Runs	39
Invalid Specimen Results.....	40
Results of Pooled Specimens (Pools of up to 24 IndividualSpecimens)	40
Negative Primary Pools	40
Positive Primary Pools — Secondary Pool Testing.....	40
Positive Primary Pool, Positive Secondary Pools — Tertiary Resolution Testing.....	41
Results of Individual Specimens.....	41
Results of Pooled Source Plasma Specimens (Pools of up to 96 Individual Specimens)	42
Negative Primary Pools	42
Positive Primary Pools — Minipool Testing.....	42
PROCEDURAL LIMITATIONS	43
PERFORMANCE CHARACTERISTICS	43
Reproducibility	43
Analytical Sensitivity — WHO HBV International Standard	46
Genotype Sensitivity and Inclusivity	47
Seroconversion Panels	49
Analytical Specificity — Potentially Cross Reactive and Interfering Microorganisms	50
Potentially Interfering Substances	51
Endogenous Interfering Substances.....	51
Exogenous Interfering Substances.....	51

CLINICAL PERFORMANCE	51
Detection of HBV DNA in Seropositive Specimens	51
Seropositive Specimens (Including 49 Acute Patients) from Commercial Vendors	52
Chronic HBV Patients.....	53
High Risk Patients.....	54
Seropositive Specimens from Seroconversion Panels	55
Seropositive Specimens from Blood Donors.....	55
Pool Reactivity in Whole Blood	56
Detection of Window Period Cases	58
Single Donation Testing Performance	59
Performance Characteristics for Source Plasma	60
Clinical Performance	60
HBV Seroconversion Panels.....	61
REFERENCES	62

INTENDED USE

The COBAS AmpliScreen HBV Test is a qualitative *in vitro* test for the direct detection of Hepatitis B Virus (HBV) DNA in human plasma.

The COBAS AmpliScreen HBV Test is intended to be used to screen donors for HBV DNA in addition to the currently recommended serology tests. This product is intended for use as a donor screening test to detect HBV DNA in plasma samples from individual human donors, including donors of Whole Blood and blood components, Source Plasma and other living donors. It is also intended for use to screen organ donors when specimens are obtained while the donor's heart is still beating. This test is not intended for use on specimens from cadaveric (non-heart-beating) donors. This test is not intended for use on samples of cord blood.

Plasma from all donors may be screened as individual samples. For donations of Whole blood and blood components, plasma samples of the donations may be tested in pools comprised of equal aliquots of not more than 24 individual samples in conjunction with licensed tests for detecting hepatitis B surface antigen (HBsAg) and antibody to hepatitis B core antigen (anti-HBc). For donations of Source Plasma, plasma samples of the donations may be tested in pools comprised of equal aliquots of not more than 96 individual samples in conjunction with licensed tests for detecting HBsAg.

This assay is not intended for use as an aid in diagnosis.

SUMMARY AND EXPLANATION OF THE TEST

Hepatitis B Virus is considered to be one of the major etiologic agents that cause chronic and acute hepatitis, cirrhosis and hepatocellular carcinoma.¹⁻⁴ HBV is one of the most infectious disease causing agents, with about 350 million chronic hepatitis B carriers worldwide.⁹ The prevalence of HBV infection in the US is approximately 200-300 thousand and in Europe it is over 1 million.^{10, 11} In the US, there has been a steady decline in hepatitis B rates with the implementation of a national vaccination strategy.²⁷ The global prevalence of chronic HBV infection, as determined by immunoserology, ranges from < 2% in western countries to ≥ 8% in Asian and African countries.¹²

HBV is a partially double-stranded circular DNA virus with a genome of approximately 3,200 bases that contains four overlapping open reading frames encoding for all viral proteins.^{1,5,6} The current risk of transfusion transmitted HBV is higher than that of HCV and HIV.^{7,8}

The presence of HBV antigens or antibodies in patients infected with HBV has led to the development of immunoserological tests that are specific for these antigens or antibodies. Implementation of these tests has reduced, but not completely eliminated, the incidence of post-transfusion hepatitis B.^{7,8} Currently, all blood units are screened for HBsAg and anti-HBc. However, it has been reported that some blood units from HBsAg and anti-HBc negative donors were associated with post-transfusion hepatitis B in the recipients^{13,14}. HBV DNA was prospectively detected in HBsAg and anti-HBc negative

blood units by Polymerase Chain Reaction (PCR).¹³ Screening of blood donations for HBV DNA will further reduce the residual transmission risk. PCR tests can also detect viremic units donated by carriers who are in the window period, early acute infection, or late resolving infection that may not be detectable by the existing immunological assays.^{13,15-17} A number of proposals have been made for performing nucleic acid tests on mini-pools composed of small aliquots from many individual units.¹⁸⁻²⁰

The COBAS AmpliScreen HBV Test uses a generic DNA specimen preparation technique in a mini-pool testing format. HBV DNA in blood donations is detected by automated amplification by PCR on the COBAS AMPLICOR™ Analyzer. The assay incorporates an Internal Control for monitoring assay performance in each individual test as well as the AmpErase® enzyme (uracil-N-glycosylase) to reduce potential contamination by previously amplified material (amplicon).

PRINCIPLES OF THE PROCEDURE

The COBAS AmpliScreen HBV Test is based on four major processes:

1. Specimen Processing
2. PCR amplification of target DNA using HBV-specific complementary primers
3. Hybridization of the amplified products to oligonucleotide probes specific to the target(s)
4. Detection of the probe-bound amplified products by colorimetric determination

Specimen Processing

Two specimen processing procedures are used with the COBAS AmpliScreen HBV Test as follows:

- Multiprep Specimen Processing Procedure for preparation of mini-pool specimens
- Standard Specimen Processing Procedure for preparation of individual specimens.

In the Standard Specimen Processing Procedure, HBV DNA is isolated directly from plasma by lysis of the virus particles with Multiprep Lysis Reagent followed by precipitation of the DNA with alcohol. In the Multiprep Specimen Processing Procedure, HBV viral particles are first pelleted from the plasma sample by high speed centrifugation, followed by lysis of the pelleted virus with a chaotropic agent (Multiprep Lysis Reagent) and precipitation of the DNA with alcohol.

The Multiprep Internal Control (MP IC), containing the HBV Internal Control, is introduced into each specimen with the Multiprep Lysis Reagent and serves as an extraction and amplification control for each processed specimen and control. The HBV Internal Control is a DNA plasmid with primer binding regions identical to those of the HBV target sequence, a randomized internal sequence of similar length and base composition as the HBV target sequence, and a unique probe binding region that differentiates the HBV Internal Control amplicon from target amplicon. These features were selected to ensure equivalent amplification of the HBV Internal Control and the HBV target DNA.

PCR Amplification

The amplification reactions are performed with the thermostable recombinant enzyme *Thermus aquaticus* DNA Polymerase (Taq pol). The reaction mixture is heated to separate double-stranded DNA. As the mixture cools, primers anneal to the target DNA and in the presence of Mg^{2+} and excess deoxynucleotide triphosphates (dNTPs), the Taq pol extends the annealed primers along the target templates to produce a double-stranded DNA molecule termed an amplicon. The COBAS AMPLICOR Analyzer automatically repeats this process for a designated number of cycles, each cycle effectively doubling the amount of amplicon DNA. The required number of cycles is preprogrammed in the COBAS AMPLICOR Analyzer.

Selective Amplification

To ensure selective amplification of nucleic acid target in the specimen and prevent amplification of pre-existing amplicon, the AmpErase[®] enzyme (uracil-N-glycosylase, UNG) is added to the COBAS AmpliScreen HBV Test. The AmpErase enzyme recognizes and catalyzes the destruction of DNA strands containing deoxyuridine²³, but not DNA containing deoxythymidine. Deoxyuridine is not present in naturally occurring DNA, but is always present in amplicon because of the use of deoxyuridine triphosphate in place of thymidine triphosphate as one of the dNTPs in the Master Mix reagent; therefore, only amplicon contain deoxyuridine. Deoxyuridine renders contaminating amplicon susceptible to destruction by AmpErase enzyme before amplification of the target DNA. AmpErase enzyme, which is included in the Master Mix reagent, catalyzes the cleavage of deoxyuridine-containing DNA rendering the DNA non-amplifiable. AmpErase is inactive at temperatures above 55°C, i.e., throughout the thermal cycling steps, and therefore does not destroy target amplicon. Following amplification, any residual enzyme is denatured by the addition of the Denaturation Solution, thereby preventing the degradation of any target amplicon.

Hybridization Reaction

Following PCR amplification, the COBAS AMPLICOR Analyzer automatically adds Denaturation Solution to the A-tubes to chemically denature the HBV amplicon and the HBV Multiprep Internal Control amplicon to form single-stranded DNA. Aliquots of denatured amplicon are then transferred to two detection cups (D-cups). A suspension of magnetic particles coated with an oligonucleotide probe specific for HBV amplicon or HBV Internal Control amplicon is added to the individual D-cups. The biotin-labeled HBV and HBV Internal Control amplicon are hybridized to the target-specific oligonucleotide probes bound to the magnetic particles. This hybridization of amplicon to the target-specific probe increases the overall specificity of the COBAS AmpliScreen HBV Test.

Detection Reaction

Following the hybridization reaction, the COBAS AMPLICOR Analyzer washes the magnetic particles in the D-cups to remove unbound material and then adds avidin-horseradish peroxidase conjugate. The avidin-horseradish peroxidase conjugate binds to the biotin-labeled amplicon. The COBAS AMPLICOR Analyzer removes unbound conjugate by washing the magnetic particles and then adds a substrate solution containing hydrogen peroxide and 3,3',5,5'-tetramethylbenzidine (TMB) to each D-cup. In the presence of hydrogen peroxide, the particle-bound horseradish peroxidase catalyzes the oxidation of TMB to form a colored complex. The absorbance is measured by the COBAS AMPLICOR Analyzer at a wavelength of 660 nm.

MATERIALS PROVIDED BY ROCHE

The COBAS AmpliScreen Multiprep Specimen Preparation and Control Kit and the COBAS AMPLICOR Wash Buffer kit are provided as stand alone kits to be used in conjunction with the COBAS AmpliScreen HBV Test, as well as the COBAS AmpliScreen HIV-1 Test, v1.5, and the COBAS AmpliScreen HCV Test.

COBAS AmpliScreen Multiprep Specimen Preparation and Control Kit

(P/N: 3302555018)

MULTIPREP/CTL

96 Test

MP (-) C

(Multiprep Negative (-) Control)

MP (+) C

(Multiprep Positive (+) Control)

MP LYS

(Multiprep Lysis Reagent)

MP DIL

(Multiprep Specimen Diluent)

MP IC

(Multiprep Internal Control)

NHP

(Negative Plasma (Human))

COBAS AmpliScreen HBV Test

(P/N: 03599779190)

HBV

96 Tests

COBAS AmpliScreen HBV Amplification Reagents

HBV AMP

HBV MMX

(HBV Master Mix)

HBV Mg²⁺

(HBV Magnesium Solution)

COBAS AmpliScreen HBV Detection Reagents

HBV DET

BH PS1

(HBV Probe Suspension 1)

BH4

(HBV Probe Suspension 2)

BI PS1

(HBV IC Probe Suspension 1)

BI4

(HBV IC Probe Suspension 2)

DN4

(Denaturation Solution)

CN4

(Avidin – Horseradish Peroxidase Conjugate)

SB3

(Substrate A)

SB

(Substrate B)

COBAS AMPLICOR Wash Buffer Kit

(P/N: 20759899123; ART: 07 5989 9; US: 83314)

WB

500 Tests

WB

(10X-Wash Concentrate)

**OTHER MATERIALS REQUIRED BUT SOLD SEPARATELY
(MAY BE PURCHASED FROM ROCHE)**

- COBAS AMPLICOR Analyzer with software version 0022B, Printer, and *Operator's Manual* for the COBAS AMPLICOR Analyzer
- COBAS AMPLICOR A-rings
- COBAS AMPLICOR D-cups
- AMPLILINK® Software, version 1.4 and *Operator's Manual* for the AMPLILINK software
- Hamilton MICROLAB® AT plus 2 Pipettor (with Hamilton SUNPLUS and RUNENDE Software, and the Roche Pooling Methods Software version 1.3), the COBAS AmpliScreen Pooling System Guide (Roche Pooling Methods Software, version 1.3 and the COBAS AmpliScreen Pooling System Guide are validated to prepare pools of equal aliquots of not more than 24 individual plasma donations using the Hamilton MICROLAB® AT plus 2 Pipettor with Hamilton SUNPLUS and RUNENDE Software).

Note: The user must validate all pooling algorithms and equipment other than those supplied by Roche.

- Sarstedt 1.5-mL tube Barcode Labels
- Hamilton Archive and Intermediate Plate Barcode Labels
- Refrigerated high speed centrifuge with fixed angle rotor (45 degrees, capacity for at least 24 x 1.5-mL tubes) with an RCF of 23,600 x g (Heraeus Centrifuge 17RS or Biofuge 28RS with HFA 22.1 rotor, Heraeus Biofuge Stratos with the 3331 rotor or equivalent)

MATERIALS REQUIRED BUT NOT PROVIDED BY ROCHE

- Microcentrifuge, (max. RCF 16,000 x g; min. RCF 12,500 x g) (Eppendorf® 5415C, HERMLE Z230M, or equivalent)
- Eppendorf® 1.25-mL Combitip® Reservoir (sterile) or equivalent
- Eppendorf Multipette® pipette or equivalent
- Ethanol, 90% or 95%, reagent grade for Molecular Biology or Histology use
- Distilled or deionized water
- Powderless disposable gloves
- Isopropyl alcohol, reagent grade
- Disposable, Sterile, Polystyrene Pipettes (5 mL, 10 mL and 25 mL)
- Sterile, RNase-free, fine-tip transfer pipettes
- Pipettors (capacity 20 µL to 1000 µL capable of providing $\pm 3\%$ accuracy and precision $\leq 5\%$) with aerosol barrier or positive displacement RNase-free tips
- Tube racks (Sarstedt P/N 93.1428 or equivalent)
- 1.5-mL sterile, non-siliconized, conical polypropylene screw-cap tubes, (Sarstedt 72.692.105 or equivalent)
- Vortex mixer
- Hamilton Slotted Deepwell Archive Plate, 2.2 mL and Sealing Capmat
- Hamilton Slotted Intermediate Plate

REAGENTS

COBAS AmpliScreen Multiprep

Specimen Preparation and Control Kit

MPSP/CTL

96 Tests

MP (–) C

8 x 0.1 mL

(Multiprep Negative (–) Control)

< 0.005% Poly rA RNA (synthetic)

EDTA

0.05% Sodium azide

MP (+) C

8 x 0.1 mL

(Multiprep Positive (+) Control)

Tris-HCl buffer

< 0.001% Non-infectious linearized plasmid DNA (microbial) containing HBV sequences

< 0.001% Non-infectious in vitro transcribed RNA (microbial) containing HCV sequences

< 0.001% Non-infectious in vitro transcribed RNA (microbial) containing HIV-1 sequences

< 0.005% Poly rA RNA (synthetic)

EDTA

0.05% Sodium azide

MP LYS

8 x 9.0 mL

(Multiprep Lysis Reagent)

Tris-HCl buffer

60% Guanidine thiocyanate

3% Dithiothreitol

< 1% Glycogen

Xn 60% (w/w) Guanidine thiocyanate

(Harmful Symbol)

MP DIL

8 x 4.8 mL

(Multiprep Specimen Diluent)

Tris-HCl buffer

< 0.005% Poly rA RNA (synthetic)

EDTA

0.05% Sodium azide

MP IC**8 x 0.1 mL****(Multiprep Internal Control)**

Tris-HCl buffer

< 0.001% Non-infectious plasmid DNA containing HBV primer binding sequences and a unique probe binding region

< 0.001% Non-infectious in vitro transcribed RNA (microbial) containing HCV primer binding sequences and a unique probe binding region

< 0.001% Non-infectious in vitro transcribed RNA (microbial) containing HIV-1 primer binding sequences and a unique probe binding region

< 0.005% Poly rA RNA (synthetic)

EDTA

<0.1% Amaranth dye

0.05% Sodium azide

NHP**16 x 1.6 mL****(Negative Plasma (Human))**

Human plasma, non-reactive by US FDA licensed tests for antibody to HIV-1/2, antibody to HCV, HBsAg, HCV RNA, HIV-1 RNA and HBV DNA.

0.1% ProClin® 300

COBAS AmpliScreen HBV Test**HBV****96 Tests****COBAS AmpliScreen HBV Amplification Reagents****HBV AMP****HBV MMX****8 x 0.7 mL****(HBV Master Mix)**

Tris HCl buffer

Glycerol

< 0.001% Taq DNA Polymerase (Taq pol, microbial)

Ammonium sulfate

0.05% dATP, dCTP, dGTP, dUTP

< 0.001% HBV-104UB and HBV-104D, HW016TBB1 primers (HBV – 104UB is biotinylated)

< 0.001% AmpErase® Uracil-N glycosylase (microbial)

< 0.5% Tween 20

0.09% Sodium azide

HBV Mg²⁺ **8 x 0.1 mL**

(HBV Magnesium Solution)

< 1% Magnesium Chloride

Amaranth dye

0.05% Sodium azide

COBAS AmpliScreen HBV Detection Reagents

HBV DET **96 Tests**

BH PS1

1 x 100 Tests

(HBV Probe Suspension 1)

MES buffer

< 0.4% Suspension of Dynabeads® (paramagnetic particles) coated with HBV-specific oligonucleotide (HBV DET) capture probe

0.09% Sodium azide

BH4

1 x 100 Tests

(HBV Probe Suspension 2)

Sodium phosphate buffer

24.9% Sodium thiocyanate

0.2% Solubilizer

(Harmful Symbol)

BI PS1

1 x 100 Tests

(HBV IC Probe Suspension 1)

MES buffer

< 0.4% Suspension of Dynabeads (paramagnetic particles) coated with HBV IC-specific oligonucleotide capture probe (SK535)

0.09% Sodium azide

BI4

1 x 100 Tests

(HBV IC Probe Suspension 2)

Sodium phosphate buffer

24.9% Sodium thiocyanate

0.2% Solubilizer

(Harmful Symbol)

DN4

1 x 100 Tests

(Denaturation Solution)

1.6% Sodium hydroxide
EDTA
Thymol blue
Xi, 1.6% (w/w) Sodium hydroxide
(Irritant Symbol)

CN4

2 x 100 Tests

(Avidin-Horseradish Peroxidase Conjugate)

Tris-HCl buffer
< 0.001% Avidin-horseradish peroxidase conjugate
Bovine serum albumin (mammalian)
Emulsit 25 (Dai-ichi Kogyo Seiyaku Co., Ltd.)
0.1% Phenol
1% ProClin 150

SB3

10 x 75 Tests

(Substrate A)

Citrate solution
0.01% Hydrogen peroxide
0.1% ProClin 150

SB

10 x 75 Tests

(Substrate B)

0.1% 3,3',5,5'-Tetramethylbenzidine (TMB)
40% Dimethylformamide (DMF)
T 40% (w/w) Dimethylformamide (DMF)
(Toxic)

R: 61-20/21-36 May cause harm to the unborn child. Harmful by inhalation and in contact with skin. Irritating to eyes.

S: 53-45 Avoid exposure – obtain special instructions before use. In case of accident or if you feel unwell, seek medical advice immediately (show the label where possible).

(Toxic Symbol)

COBAS AMPLICOR Wash Buffer

WB

500 Tests

WB

(10X-Wash Concentrate)

2 x 250 Tests

< 2% Phosphate buffer

< 9.0% Sodium chloride

EDTA

< 2% Detergent

0.5% ProClin 300

STORAGE INSTRUCTIONS

- A. Room temperature is defined as 15-30°C.
- B. Do not freeze reagents.
- C. Store the following reagents at 2-8°C. Unopened, these reagents are stable until the expiration date indicated.

MP LYS, MP IC, MP (+) C, MP (-) C, MP DIL, and NHP

HBV MMX, HBV Mg²⁺

BH PS1, BH4, BI PS1, BI4

CN4, SB3, SB

- D. Store **DN4** and **WB** at 2-25°C. **DN4** and **WB** are stable until the expiration dates indicated.
- E. Do not expose **SB3**, **SB** or Working Substrate to metals, oxidizing agents, or direct sunlight.
- F. The following reagents are one-time use. Discard any unused portion.

MP IC, MP (+) C, MP (-) C, MP DIL, and NHP

HBV Mg²⁺, and SB

PRECAUTIONS

For *In Vitro* Diagnostic Use.

- A. Specimens may be infectious. Use Universal Precautions when performing the assay.²⁴⁻²⁵ Only personnel proficient in the use of the COBAS AmpliScreen System and trained in handling infectious materials should perform this procedure. Thoroughly clean and disinfect all work surfaces with a freshly prepared solution of 0.5% sodium hypochlorite in distilled or deionized water. Follow by wiping down the surface with 70% ethanol.
- B. **CAUTION: The Negative Human Plasma of this kit contains human blood products non-reactive by US FDA licensed tests for antibody to HIV-1/2, antibody to HCV, HBsAg, HCV RNA, HIV-1 RNA and HBV DNA. No known test method can offer complete assurance that products derived from human blood will not transmit infectious agents. All human blood-sourced materials should be considered potentially infectious and should be handled with Universal Precautions.** If spillage occurs, immediately disinfect, then wipe up with a 0.5% (final concentration) sodium hypochlorite solution (diluted bleach) or follow appropriate site procedures.
- C. Use routine laboratory precautions. Do not pipette by mouth. Do not eat, drink or smoke in designated work areas. Wear disposable gloves, laboratory coats and eye protection when handling specimens and kit reagents. Wash hands thoroughly after handling specimens and kit reagents.
- D. This product contains sodium azide as a preservative. Do not use metal tubing for reagent transfer. If solutions containing azide compounds are disposed of in a plumbing system, they should be diluted and flushed with generous amounts of running water. These precautions are recommended to avoid accumulation of deposits in metal piping in which explosive conditions could develop.

-
- E. ***Heparin has been shown to inhibit PCR. Do not use heparinized plasma with this procedure.***
- F. Use only supplied or specified required disposables to ensure optimal assay performance.
- G. Screw-cap tubes must be used for specimen and control preparation to prevent splashing and potential cross-contamination of specimens and controls. ***Do not use snap cap tubes.***
- H. Adequately vortex, where specified, to ensure optimal assay performance.
- I. Handle all materials containing specimens or controls according to Good Laboratory Practices in order to prevent cross-contamination of specimens or controls.
- J. Before use, visually inspect each reagent bottle to ensure that there are no signs of leakage and/or abnormal color. If there is any evidence of leakage and/or abnormal color, do not use that bottle for testing.
- K. Dispose of all materials that have come in contact with specimens and reagents in accordance with country, federal, state and local regulations.
- L. Do not use a kit after its expiration date. DO NOT interchange, mix, or combine reagents from kits with different master lot numbers. Do not use expired reagents.
- M. Material Safety Data Sheets (MSDS) are available on request.
- N. Supplies and equipment must be dedicated to each pre-amplification activity and should not be used for other activities or moved between areas. **Fresh, clean gloves must be worn in each area and must be changed before leaving that area.**
Equipment and supplies used for reagent preparation must not be used for specimen preparation activities or for pipetting or processing amplified DNA or other sources of target DNA. Post-amplification supplies and equipment must remain in the Post-Amplification Area at all times.
- O. Avoid contact of **MP LYS, HBV MMX, HBV Mg²⁺, BH4, BI4, DN4, CN4, SB3, SB** and Working Substrate (mixed **SB3** and **SB** reagent) with the skin, eyes or
-

mucous membranes. **If contact does occur, immediately wash with large amounts of water, otherwise, burns can occur.** If these reagents are spilled, dilute with water before wiping dry. **Do not allow MP LYS, which contains guanidine thiocyanate, or BH4 and BI4, which contain sodium thiocyanate, to contact sodium hypochlorite (bleach) solution. This mixture can produce a highly toxic gas.**

- P. SB and Working Substrate contain dimethylformamide, which has been reported to be toxic in high oral doses, and may be harmful to the unborn child. Skin contact, inhalation of fumes and ingestion should be avoided. If skin contact occurs, wash thoroughly with soap and water and seek medical advice immediately.**
- Q.** Refer to “Precautions” in other COBAS AmpliScreen package inserts, COBAS Pooling System Guide, and the Operator’s Manuals for the AMPLILINK and COBAS AMPLICOR Analyzer.
- R.** Closely follow procedures and guidelines provided to ensure that the specimen and control preparation is performed correctly. Any deviation from the given procedures and guidelines may affect optimal assay performance.

REAGENT PREPARATION

A. MP IC, MP (+) C, MP (–) C, MP DIL and NHP

1. Warm MP IC, MP(+) C, MP(-) C, MP DIL and NHP to room temperature before use by using a 37°C incubator or on laboratory bench top.

B. Working Lysis Reagent

1. Warm MP LYS to 25 – 37°C to dissolve precipitate (maximum 30 minutes). Mix thoroughly until the crystals are dissolved. Prior to use, visually verify that crystals are dissolved and examine each bottle of MP LYS against a white background for the appearance of a yellow color or leakage. If there is any

yellow color or signs of leakage do not use that bottle for testing. Contact your local Roche office for replacement.

2. Vortex MP IC briefly before use. Tap vial to collect the solution in the base. Pipette 100 μ L MP IC into 1 bottle MP LYS. Cap the MP LYS bottle and vortex briefly. The pink color confirms that the MP IC has been added to the MP LYS. Discard the remaining MP IC.
3. Store Working Lysis Reagent at room temperature. Use within 4 hours of preparation.

C. Working Amplification Master Mix

1. Prepare Working Master Mix in a template-free area (e.g., in a dead air box).
Reagent preparation area must be clean and disinfected in accordance with methods outlined in Precautions (Item A). Failure to do so may result in reagent contamination.
2. Pipette 100 μ L HBV Mg²⁺ into 1 bottle HBV MMX. Recap HBV MMX bottle and mix well by inverting 10 – 15 times. The pink color confirms that the HBV Mg²⁺ has been added to the HBV MMX. Discard the remaining HBV Mg²⁺. Do not vortex the Working Master Mix. These reagents do not need to be at room temperature before use.
3. Store at 2 – 8°C and use within 4 hours of preparation.

D. Working Probe Suspension Detection Reagents

1. Prepare Working HBV Probe Suspension: Mix BH PS1 well by vortexing briefly to suspend the microparticles. Pipette 2.5 mL BH PS1 into one BH4 cassette.
2. Prepare Working IC Probe Suspension: Mix BI PS1 well by vortexing briefly to suspend the microparticles. Pipette 2.5 mL BI PS1 into one BI4 cassette.
3. Both Working Probe Suspension Detection Reagents are stable for 30 days at 2 – 8°C. Working Reagents can be used for a maximum of six instrument

cycles (12 hours per cycle). Mixing occurs automatically on the COBAS AMPLICOR Analyzer.

4. Store Working Probe Suspension Detection Reagents at 2 – 8°C between instrument cycles. Remove from refrigerator 30 minutes before use on the COBAS AMPLICOR Analyzer.

E. DN4 — Denaturation Reagent and CN4 Conjugate Reagent

1. Once opened, DN4 and CN4 are stable for 30 days at 2 – 8°C, or until the expiration date, whichever comes first. Both DN4 and CN4 can be used for a maximum of six instrument cycles (12 hours per cycle).
2. Store DN4 and CN4 at 2 – 8°C between instrument cycles. Remove from refrigerator 30 minutes before use on the COBAS AMPLICOR Analyzer.

F. Working Substrate Reagent

1. Working Substrate must be prepared each day by pipetting 5 mL **SB** into one **SB3** cassette. Pipette up and down at least 5 times to mix.
2. Working Substrate is stable on the COBAS AMPLICOR Analyzer for a maximum of 16 hours.
3. Do not expose **SB3**, **SB**, or Working Substrate to metals, oxidizing agents or direct light.

G. Wash Buffer Reagent

1. Examine WB before dilution and if necessary, warm at 30 – 37°C to dissolve any precipitate. Add 1 volume of WB to 9 volumes of distilled or deionized water. Mix well. Keep a minimum of 3 – 4 liters of Working Wash Buffer (1X) in the Wash Buffer Reservoir of the COBAS AMPLICOR Analyzer at all times.

2. Working Wash Buffer (1X) should be stored at 2 – 25°C in the COBAS AMPLICOR Wash Buffer Reservoir and is stable for 2 weeks from the date of preparation.


H. 70% Ethanol

1. Prepare 70% ethanol fresh daily.
2. One mL 70% ethanol is needed for each specimen and control processed. For example, mix 11.7 mL 90% ethanol and 3.3 mL of distilled or deionized water for every 12 specimens and controls to be processed.

SPECIMEN COLLECTION, STORAGE AND POOLING

Note: Handle all specimens as if they are potentially infectious agents.

- A. EDTA, CPD, CPDA-1, CP2D, ACD-A and 4% Sodium Citrate may be used with the COBAS AmpliScreen HBV Test. Follow sample tube manufacturer's instructions.
- B. Blood collected in EDTA may be stored at 2 – 30°C for up to 72 hours from time of draw, followed by an additional two days at 2 – 8°C. For storage longer than five days, remove the plasma from the red blood cells by centrifugation at 800 – 1600 x g for 20 minutes. Following removal, plasma may be stored at 2 – 8°C for an additional seven days. Alternatively, plasma may be stored at $\leq -18^{\circ}\text{C}$ for up to one month.

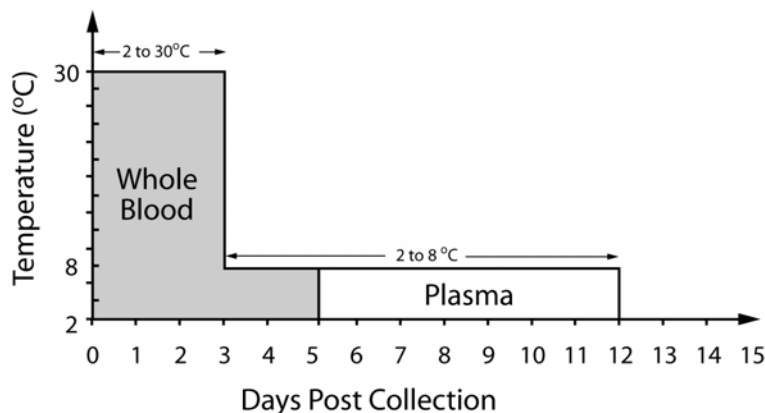

- C. Blood collected in CPD, CPDA-1, or CP2D may be stored for up to 72 hours at 1 – 24°C. Following centrifugation of the CPD, CPDA-1, or CP2D samples at 800-1600 x g for 20 minutes, plasma may be stored at 1 – 6°C for an additional 7 days from the date the plasma was removed from the red blood cells. Plasma separated from the cells may be stored at $\leq -18^{\circ}\text{C}$ for up to one month.
- D. ACD-A or 4% sodium citrate anticoagulated apheresis plasma can be stored at 1 – 6°C for up to 6 hours, followed by subsequent storage at $\leq -18^{\circ}\text{C}$ for up to one month.
- E. **Do not freeze whole blood.**
- F. **Heparin has been shown to inhibit PCR. Use of heparinized specimens is not recommended.**
- G. Warm pooled or individual specimens to room temperature before using.
- H. Covered Archive Plates may be stored at 2 – 8°C for up to 7 days from the date the plasma was removed from the red blood cells.
- I. No adverse effect on assay performance was observed when plasma specimens were subjected to three freeze-thaw cycles.
- J. Thaw frozen specimens at room temperature before using.
- K. The user should validate other collection and storage conditions. If specimens are to be shipped, they should be packaged and labeled in compliance with applicable federal and international regulations covering the transport of clinical specimens and etiologic agents.²⁶
- L. **False positive results may occur if cross contamination of specimens is not adequately controlled during specimen handling and processing.**
- M. **SPECIMEN POOLING**
 - 1. The COBAS AmpliScreen Pooling System performs barcode scanning and pooling operations that combine aliquots from 24 individual specimens into a

single Primary Pool that is used for testing. The pooling algorithm requires preparation of Secondary Pools as well as individual specimens for follow-up testing in the event a Primary Pool tests positive. If less than 24 specimens are available, testing is performed using the individual specimens.

2. For Source Plasma, the Hamilton performs barcode scanning and pooling operations that combine aliquots from 96 individual specimens into a single Primary Pool that is used for testing. Positive Primary Pools are traced to the positive individual using an overlapping pool testing matrix. Minipools are prepared from the eight individual specimens for columns 1 – 12 and from the 12 individual specimens for rows 1 – 8. The positive specimen is identified by the intersection of the positive column and positive row. Confirmatory testing is conducted on the implicated specimen using Standard Specimen Processing Procedure. (Hamilton MICROLAB AT plus 2 Pipettor with SUNRISE PLUS v3.3 software was used to prepare pools of up to 96 equal aliquots of plasma during the clinical trials).

Note: The user must validate all pooling algorithms and equipment other than those supplied by Roche.

PROCEDURAL NOTES

A. Run Size

1. Each kit contains reagents sufficient for eight 12-specimen runs, which may be performed separately or simultaneously. At least one preparation of the COBAS AmpliScreen Multiprep Negative (–) Control and one preparation of the COBAS AmpliScreen Multiprep Positive (+) Control must be included in each run (see “Quality Control” section).
2. The Specimen Preparation and Amplification Reagents are packaged in eight single-use bottles. The Multiprep Negative (–) and Multiprep Positive (+) Controls are packaged in single-use vials. For the most efficient use of

reagents, specimens and controls should be processed in batches that are multiples of 12.

3. The use of sterile gauze, when uncapping specimen tubes may reduce the potential for cross contamination between specimens.

B. Equipment

1. Prepare the COBAS AMPLICOR Analyzer and AMPLILINK Data Station for use according to instructions in the Operator's Manual for the AMPLILINK software and the Operator's Manual for the COBAS AMPLICOR Analyzer.
2. Prepare the Hamilton Microlab AT plus 2 System and SUNPLUS Data Station for use according to instructions in the Operator's Manuals.
3. Pre-cool the high-speed centrifuge and rotor to 2 – 8°C. See operating instructions for the high-speed centrifuge for details.
4. Perform manufacturer recommended maintenance and calibration on all instruments, including pipettors, to ensure proper functioning.

C. Reagents

1. All reagents, except HBV MMX and HBV Mg²⁺ must be at room temperature before use. Visually examine reagents for sufficient volume before beginning the test procedure. See section "Reagent Preparation" for specific reagent storage conditions.
2. Add all reagents using a pipettor capable of delivering specified volume with $\pm 3\%$ accuracy and a precision of $\leq 5\%$ CV. Check pipettor functionality and calibrate as recommended by pipettor manufacturer.
3. Prepare Working Master Mix in a template-free area (e.g., in a dead air box). Reagent preparation area must be clean and disinfected in accordance with methods outlined in "Precautions" (Item A). Failure to do so may result in reagent contamination.

4. Prepare 70% ethanol fresh each day.
5. Check expiration date of opened or Working Reagents before loading on the COBAS AMPLICOR Analyzer.
6. Check to ensure that all reagents used are of the same master lot of kit reagents.

D. Workflow

1. To minimize the possibility of laboratory areas becoming contaminated with amplicon, the laboratory area should be separated into several distinct areas organized around Pre-Amplification and Post-Amplification. Personnel should use proper anti-contamination safeguards when moving between areas.
2. The Pre-Amplification Area should have a template-free area for preparation of Working Master Mix and an amplicon free area for specimen and control preparation.
3. The Post-Amplification Area should have a COBAS AMPLICOR Analyzer(s) and AMPLILINK Data Station(s) with additional area for preparing Working Amplification and Detection Reagents.
4. Pipettors and other supplies should be dedicated to a specific area. Specimens, equipment and reagents should not be returned to the area where a previous step was performed.

E. Temperature

1. Room temperature is defined as 15° to 30°C.

F. Vortexing

1. Proper vortexing during specimen preparation is important to ensure homogeneous mixture after additions of reagents.

G. Pipetting

1. Pooled or individual plasma specimens must be at room temperature before pipetting.
2. Use a clean pipette tip or disposable transfer pipette with each specimen or control. Use aerosol barrier or positive displacement RNase-free tips.
3. Confirm that all pipettors are correctly set to dispense the specified volumes in accordance with the specimen preparation procedures and guidelines.

H. Specimen Processing

1. Screw-cap tubes must be used for specimen and control preparation to prevent splashing and potential cross-contamination of specimens and controls. Do not use snap cap tubes.
2. Avoid contaminating gloves when manipulating specimens.
3. Specimens and controls should be prepared in a laminar flow hood. *Failure to do so may result in specimen contamination.* Specimen and control preparation area must be cleaned and disinfected in accordance with methods outlined in “Precautions” (Item A).

I. Decontamination

1. Thoroughly clean and disinfect all work surfaces with a freshly prepared solution of 5% sodium hypochlorite in distilled or deionized water. Follow by wiping down the surface with 70% ethanol.

INSTRUCTIONS FOR USE

The Multiprep Specimen Processing Procedure is used for extracting nucleic acid from pooled specimens. The Standard Specimen Processing Procedure is used for extracting nucleic acid from individual specimens.

The Multiprep and the Standard Specimen Processing Procedures are generic nucleic acid extraction procedures and can be used for the extraction of HIV-1 RNA, HCV RNA, and/or HBV DNA. A single extraction is sufficient for multiple assays. Workflow can be performed on the same day or over multiple days under the following conditions:

Amplification, Hybridization and Detection of Stored Processed Specimens

Amplification, hybridization and detection can occur on the same day as specimen processing or on a subsequent day. If amplification, hybridization and detection are to be done on a subsequent day, perform the Multiprep Specimen Processing Procedure described in steps B1 through B21 or the Standard Specimen Processing Procedure described in steps B22 through B38. Store the processed specimens and controls as indicated. On the subsequent day, begin with Step A (Reagent Preparation — Working Master Mix), thaw processed specimens and controls at room temperature, and continue with Step B39.

Hybridization and Detection of Stored Denatured Amplicon

Hybridization and detection of the denatured amplicon may occur on the same day as amplification or on a subsequent day. If hybridization and detection are to be done on a subsequent day, the denatured amplicon may be left on-board the COBAS AMPLICOR Analyzer for not more than 24 hours before starting the hybridization and detection steps. Alternatively, the denatured amplicon may be stored at 2-8°C for not more than five days before starting the hybridization and detection steps.

A. Reagent Preparation — Working Master Mix

Performed in: Pre-Amplification — Reagent Preparation Area (e.g., dead air box)

- A1. Determine the appropriate number of A-ring(s) needed for specimen and control testing.
- A2. Place the A-ring(s) on the A-ring holder(s).
- A3. For each A-ring, prepare one Working Master Mix.
- A4. Pipette 50 μ L Working Master Mix into each A-tube. Discard unused Working Master Mix. Do not close the covers of the A-tubes at this time.
- A5. Place the A-ring containing Working Master Mix in a sealable bag and seal the plastic bag. Record the assay name (HBV) and the time the Working Master Mix was prepared.
- A6. Store the A-ring(s) containing Working Master Mix at 2-8°C until specimen and control preparation is completed. The A-rings with Working Master Mix must be used within 4 hours of preparation.
- A7. Decontaminate area. See “*Procedural Notes*”, Item I.

B. Specimen and Control Preparation

Performed in Pre-Amplification — Specimen and Control Preparation Area

Multiprep Specimen Processing Procedure

- B1. Pipette 1000 μ L of each pool into an appropriately labeled screw-cap tube using the COBAS AmpliScreen Pooling System, a hand-held pipettor or other user-validated method. Cap the tubes.
- B2. Vortex NHP briefly.
- B3. For each Negative and Positive Control pipette 1000 μ L NHP into an appropriately labeled screw-cap tube. Cap the tubes.
- B4. Use a permanent marker to make an orientation mark on each tube.

-
- B5. Place the specimen and control tubes into the pre-cooled high-speed centrifuge with the orientation marks facing outward, so that the orientation marks will align with the pellets formed during centrifugation.
- B6. Centrifuge specimens and control tubes at 23,000 - 24,000 x g for 60 ± 4 minutes at 2-8°C. The pellet will form on the outer wall as indicated by the orientation mark.

Note: The 60 ± 4 minutes begins when the centrifuge reaches 23,000 - 24,000 x g.

- B7. Remove the tubes from the centrifuge and remove the caps. Slowly aspirate 900 µL of the supernatant from each centrifuged tube leaving approximately 100 µL of supernatant. Avoid contact with the pellet. Discard the supernatant and pipette tip appropriately. Use a fresh pipette tip for each tube.
- B8. Prepare a Working Lysis Reagent bottle for every batch of 12 specimens and controls to be processed.
- B9. Pipette 600 µL Working Lysis Reagent into each specimen and control tube. Cap and vortex tubes briefly.
- B10. Prepare Controls as follows:
- a. Negative Control

Vortex MP (–) C briefly. Tap vial to collect the solution in the base.
Pipette 20 µL MP (–) C to the tube labeled “MP (–) C” containing Working Lysis Reagent and NHP. Cap the tube and vortex briefly.
 - b. Positive Control

Vortex MP (+) C briefly. Tap vial to collect the solution in the base.
Pipette 20 µL MP (+) C to the tube labeled “MP (+) C” containing Working Lysis Reagent and NHP. Cap the tube and vortex briefly.
- B11. Incubate all tubes for 10 to 15 minutes at room temperature after adding Working Lysis Reagent to the last tube. After the incubation period, briefly vortex all tubes.

- B12. Pipette 700 μ L of isopropanol into each tube. Cap the tubes and vortex briefly.
- B13. Place the tubes into a microcentrifuge with the orientation marks facing outward to align with the pellets that will form. Centrifuge at $14,250 \pm 1750 \times g$ for 15-20 minutes at room temperature.
- B14. Slowly aspirate the supernatant from each tube. Remove as much liquid as possible without disturbing the pellet.
- B15. Pipette 1.0 mL of 70% ethanol into each tube. Cap the tubes and vortex briefly.
- B16. Place the tubes into a microcentrifuge with the orientation marks facing outward to align with the pellets that will form. Centrifuge at $14,250 \pm 1750 \times g$ for 5-10 minutes at room temperature.
- B17. Slowly aspirate the supernatant from each tube using a fine-tip disposable transfer pipette. Remove as much liquid as possible without disturbing the pellet. Use a new transfer pipette for each tube.
- B18. Using a new transfer pipette for each tube, repeat Step B17 to remove as much of the remaining supernatant as possible without disturbing the pellet. ***Residual ethanol can inhibit amplification.***
- B19. Pipette 200 μ L MP DIL into each tube. Use a pipette tip to break apart the pellet. This can be done by aspirating 30-40 μ L of the diluent in the tip and scraping the sides and base of the tube in an up/down motion for at least 10 seconds and dispensing 30-40 μ L. Cap the tubes and vortex briefly to resuspend the extracted DNA. Note that some insoluble material may remain.
- B20. At this point amplification of the processed specimens and controls must be started within 2 hours. If not, the processed specimens and controls can be stored at -70°C or colder for up to one month. ***Thawing should be completed within one hour at room temperature.***
- B21. Proceed to step B39, Loading the A-ring.
-

Standard Specimen Processing Procedure

- B22. Pipette 200 μ L of each specimen into an appropriately labeled screw-cap tube using the COBAS AmpliScreen Pooling System, a hand-held pipettor or other user-validated method. Cap the tubes.
- B23. Vortex NHP briefly.
- B24. For each Negative and Positive Control pipette 200 μ L NHP into appropriately labeled screw-cap tubes. Cap the tubes.
- B25. Use a permanent marker to make an orientation mark on each tube.
- B26. Prepare a Working Lysis Reagent bottle for every 12 specimens and controls to be processed.
- B27. Pipette 600 μ L Working Lysis Reagent into each tube. Cap and vortex tubes briefly.
- B28. Prepare Controls as follows:
- a. Negative Control

Vortex MP (–) C briefly. Tap vial to collect the solution in the base.
Pipette 20 μ L MP (–) C into the tube labeled “MP (–) C” containing Working Lysis Reagent and NHP. Cap the tube and vortex briefly.
 - b. Positive Control

Vortex MP (+) C briefly. Tap vial to collect the solution in the base.
Pipette 20 μ L MP (+) C into the tube labeled “MP (+) C” containing Working Lysis Reagent and NHP. Cap the tube and vortex briefly.
- B29. Incubate all tubes for 10-15 minutes at room temperature after adding Working Lysis Reagent to the last tube. After the incubation period, briefly vortex all tubes.
- B30. Pipette 800 μ L of isopropanol into each tube. Cap the tubes and vortex briefly.

-
- B31. Place the tubes into a microcentrifuge with the orientation marks facing outward to align with the pellets that will form. Centrifuge at $14,250 \pm 1750 \times g$ for 15-20 minutes at room temperature.
- B32. Slowly aspirate the supernatant from each tube. Remove as much liquid as possible without disturbing the pellet.
- B33. Pipette 1.0 mL of 70% ethanol into each tube. Cap the tubes and vortex briefly.
- B34. Place the tubes into a microcentrifuge with the orientation marks facing outward to align with the pellets that will form. Centrifuge at $14,250 \pm 1750 \times g$ for 5-10 minutes at room temperature.
- B35. Slowly aspirate the supernatant from each tube using a fine-tip disposable transfer pipette. Remove as much liquid as possible without disturbing the pellet. Use a new transfer pipette for each tube.
- B36. Using a new transfer pipette for each tube, repeat Step B35 to remove as much of the remaining supernatant as possible without disturbing the pellet. ***Residual ethanol can inhibit amplification.***
- B37. Pipette 200 μ L MP DIL into each tube. Use a pipette tip to break apart the pellet. This can be done by aspirating 30-40 μ L of the diluent in the tip and scraping the sides and base of the tube in an up/down motion for at least 10 seconds and dispensing 30-40 μ L. Cap the tubes and vortex briefly to resuspend the extracted DNA. Note that some insoluble material may remain.
- B38. At this point amplification of the processed specimens and controls must be started within 2 hours. If not, the processed specimens and controls can be stored at -70°C or colder for up to one month. ***Thawing should be completed within one hour at room temperature.***

Loading the A-ring

- B39. Create an A-ring worklist record for each A-ring to identify the A-tube with the appropriate control or specimen to be pipetted.

- B40. If processed specimens and controls were stored frozen, thaw at room temperature before proceeding. Briefly vortex the processed specimens and controls.
- B41. Pipette 50 µL of each processed specimen and control into the appropriate A-tube containing HBV Working Master Mix. Immediately cap the A-tube and repeat this step for all the 12 A-tubes to complete the A-ring loading. Use the A-ring worklist record to ensure the appropriate specimen or control is added to the correct A-tube position for each A-ring.
- B42. Transfer the A-ring with sealed tubes containing the processed specimens and controls in Working Master Mix to the Amplification/Detection Area. Proceed to Part C.

Note: Amplification must begin within 45 minutes from when the first specimen or control in the A-ring is added to the Working Master Mix.

C. Amplification and Detection

Performed in Post-Amplification — Amplification/Detection Area

- C1. Perform Daily Instrument Maintenance as outlined in the Operator's Manual for the COBAS AMPLICOR Analyzer including:
- a. Wipe D-cup handler tip with a lint-free moist cloth and dry.
 - b. Wipe initialization post with a lint-free moist cloth and dry.
- C2. Before each run:
- a. Check waste container and empty if necessary.
 - b. Check Wash Buffer Reservoir and add prepared Wash Buffer if necessary.
 - c. Replace used D-cup racks.
 - d. Prime the COBAS AMPLICOR Analyzer.

C3. Instrument Loading and System Operation

- a. Prepare enough of the following detection reagent cassettes to complete the workload: Working HBV Probe Suspension Reagent (BH4), Working IC Probe Suspension Reagent (BI4), Working Substrate (SB3), Denaturation Reagent (DN4), and Conjugate Reagent (CN4).
- b. Place the BH4 and BI4 cassettes in the test-specific reagent rack.
- c. Place DN4, CN4 and SB3 cassettes in the generic reagent rack. Record on the cassette the date when each cassette was opened.
- d. Identify the reagent racks as generic or test specific using the COBAS AMPLICOR Analyzer barcode scanner for the AMPLILINK software, as described in the Operator's Manual for AMPLILINK software.
- e. Configure the reagent racks by entering the reagent positions and lots using the COBAS AMPLICOR Analyzer barcode scanner for the AMPLILINK software, as described in the Operator's Manual for AMPLILINK software.
- f. Load the reagent racks onto the analyzer using the COBAS AMPLICOR Analyzer barcode scanner for the AMPLILINK software, as described in the Operator's Manual for AMPLILINK software. Make sure that each reagent cassette is in its assigned position and that each cassette fits tightly into its rack.
- g. Place the D-cup rack on the D-cup platform. Two D-cups are required for each A-tube and two D-cups are required for each Working Substrate cassette to allow for blanking by the COBAS AMPLICOR Analyzer, as described in the Operator's Manual for the COBAS AMPLICOR Analyzer.
- h. Place the A-ring into the thermal cycler segment of the COBAS AMPLICOR Analyzer and close the cover on the thermal cycler segment.

- i. Load the A-ring into the COBAS AMPLICOR Analyzer using the Analyzer barcode scanner for the AMPLILINK software, as described in the Operator's Manual for AMPLILINK software.
- j. Create an A-ring order, using the AMPLILINK software, as described in the Operator's Manual for AMPLILINK software. Use the A-ring worklist record created for specimen processing to assist in entering the A-ring order.
- k. Repeat steps h. through j. above to load a second A-ring on the COBAS AMPLICOR Analyzer.
- l. Start the COBAS AMPLICOR Analyzer as described in the Operator's Manual for AMPLILINK software.
- m. Wait for the COBAS AMPLICOR Analyzer to indicate that the load check has passed.

Note: The required quantity of each detection reagent is automatically calculated by the COBAS AMPLICOR Analyzer during the Load Check to determine if sufficient reagents are available for the requested tests.

- n. The COBAS AMPLICOR Analyzer automatically performs amplification and detection. Results are expressed as absorbance values at 660 nm and as positive or negative.
- o. As a Quality Control measure, the AMPLILINK A-ring Results Report and the Run Log may be printed (e.g. daily, weekly or monthly) and retained along with the respective A-ring worklist. A selection of A-ring worklist records should be periodically compared with the AMPLILINK A-ring Results Report to verify that the A-ring ID, instrument serial number, and specimen IDs are identical. Reconcile the Run Log with the selected A-ring worklist to account for all A-ring IDs associated with the run. If there are discrepancies, perform follow-up investigation.

QUALITY CONTROL PROCEDURES

1. At least one Multiprep (–) Control and one Multiprep (+) Control must be processed with each A-ring.

- a. Negative Control

The absorbance for the MP (–) C should be less than 0.2 at 660 nm and its associated MP IC should be greater than or equal to 0.2 for the Negative Control to be valid. If the absorbance value for the MP (–) C is greater than or equal to 0.2 and/or its associated MP IC is less than 0.2, the entire A-ring is invalid, and the entire test procedure for that A-ring (specimen and control preparation, amplification and detection) must be repeated.

- b. Positive Control

The absorbance for the MP (+) C should be greater than or equal to 1.0 at 660 nm and its associated MP IC should be greater than or equal to 0.2 at 660 nm for the Positive Control to be valid. If the absorbance value for the MP (+) C is less than 1.0 and/or its associated MP IC is less than 0.2, the entire A-ring is invalid, and the entire test procedure for that A-ring (specimen and control preparation, amplification and detection) must be repeated.

Summary of Control Acceptance Criteria

	HBV Result		IC Result	
	A ₆₆₀	Comment	A ₆₆₀	Comment
Negative Control	< 0.2	Negative	≥ 0.2	Valid
Positive Control	≥ 1.0	Positive	≥ 0.2	Valid

2. Flags and comments may be generated by the COBAS AMPLICOR Analyzer during a run. The Operator must check the run printout(s) for flags and comments to verify that the run is valid. Refer to the *Operator's Manual* for the

AMPLILINK software and the *Operator's Manual* for the COBAS AMPLICOR Analyzer for interpretation of flags and comments.

a. External Control

If an External Control (i.e., an additional run control other than the Multiprep (+) Control or Multiprep (–) Control) is required by the laboratory, the External Control should meet regulatory requirements for such controls. The absorbance of the HBV External Control should be equal to or greater than 0.2 at 660 nm, irrespective of the MP IC absorbance. If the absorbance of the HBV External Control does not meet the above criterion, the negative results for specimens in the associated run may be invalidated. However, positive results for specimens in such a run should not be invalidated solely on the basis of the results obtained for an External Control; those positive results should remain the test of record. The laboratory should follow its established Standard Operating Procedure for the appropriate action.

INTERPRETATION OF RESULTS

1. Flags and comments may be generated by the COBAS AMPLICOR Analyzer during a run. The Operator must check the run printout(s) for flags and comments to verify that the run is valid. Refer to the *Operator's Manual* for the AMPLILINK software and the *Operator's Manual* for the COBAS AMPLICOR Analyzer for interpretation of flags and comments.

2. Specimen Results

Two absorbance values are obtained for each specimen: one for the HBV target and one for the internal control MP IC. For a specimen with an absorbance less than 0.2, the MP IC absorbance for that specimen must be greater than or equal to 0.2 at 660 nm for a valid negative specimen test result. If the absorbance for the

HBV target is greater than or equal to 0.2, the MP IC result is disregarded and the test result is valid and positive.

3. For a valid run, results are interpreted as follows:

HBV Result		IC Result		Interpretation
A ₆₆₀	Comment	A ₆₆₀	Comment	
< 0.2	NEGATIVE	≥ 0.2	VALID	Specimen is negative for HBV DNA.
< 0.2	NEGATIVE	< 0.2	INVALID	Invalid result. Repeat entire test procedure for invalid specimen.
≥ 0.2	POSITIVE	ANY	VALID	Specimen is positive for HBV DNA.

Invalid Test Runs

When invalid Positive or Negative Control results are obtained on an A-ring, that A-ring is invalid. Repeat the entire test procedure for the associated specimens (including specimen and control preparation, amplification and detection) in the A-ring by processing another aliquot of the original plasma specimens.

With the exception of instrument failures subsequent to denaturation of amplicon, an instrument failure during a test run, as indicated by system error messages, also constitutes an invalid test run. In such instances, repeat the test procedure for the associated controls and specimens (amplification and detection) in the run by processing another aliquot of the processed specimen.

For instrument failures subsequent to successful denaturation of amplicon, it is not necessary to repeat the entire test procedure for the associated specimens. In such instances, the denatured amplicon may be redetected by the COBAS AMPLICOR Analyzer. The denatured amplicon may be left on the COBAS AMPLICOR Analyzer for not more than 24 hours before continuing with the hybridization and detection steps. Alternatively, the denatured amplicon may be stored at 2 – 8°C for not more than five days before continuing with the hybridization and detection steps.

Invalid Specimen Results

For specimen(s) that are invalid, repeat entire test procedure in single on the remaining replicate tubes. The test result for the pool or individual specimen is based only on the repeat valid test result. If the last available replicate of a pooled specimen gives an invalid result, each individual specimen in that pool should be tested. If an individual specimen gives an invalid result, the test result for that individual specimen should be considered invalid for HBV DNA.

Results of Pooled Specimens (Pools of up to 24 Individual Specimens)

The testing algorithm for testing of pooled specimens for the COBAS AmpliScreen HBV Test requires a single level of testing for Primary Pools that are negative for HBV DNA and three levels of testing (Primary Pool, Secondary Pool and tertiary resolution) for Primary Pools that are positive for HBV DNA.

Negative Primary Pools

When the Primary Pool is negative, report the results for all associated individual specimens in that Primary Pool as “HBV DNA Negative”.

Positive Primary Pools — Secondary Pool Testing

When the Primary pool is positive, each of the 4 Secondary Pools prepared from the Archive Plate using specimens comprising the Primary Pool is tested by the COBAS AmpliScreen HBV Test, using the Multiprep Specimen Processing Procedure.

- If one or more of the Secondary Pools tests positive, report the results for the specimens in the negative Secondary Pools as “HBV DNA negative”. For positive Secondary Pools, proceed to the section entitled “Positive Primary Pool, Positive Secondary Pools — Tertiary Resolution Testing.”
- If all 4 Secondary Pools are negative, the individual specimens in that Primary Pool may be reported as “HBV DNA Negative.”

- As part of an overall Quality Assurance program, you may wish to conduct additional testing to determine the cause of the initial positivity of the Primary Pool.

Positive Primary Pool, Positive Secondary Pools — Tertiary Resolution Testing

For a positive Secondary Pool, test each of the individual specimens in that Secondary Pool. The individual specimens must be processed using the Standard Specimen Processing procedure.

- If one or more of the individual specimens is positive, the positive specimen(s) is reported as “HBV DNA positive” and the remaining negative specimens associated with the positive Secondary Pool are reported as “HBV DNA negative”.
- If all of the individual specimens in that Secondary Pool test negative, the specimens in the Secondary Pool may be reported as “HBV DNA negative”.
- As part of an overall Quality Assurance program, you may wish to conduct additional testing to determine the cause of the initial positivity of the Primary and Secondary Pools.

Results of Individual Specimens

- If an individual specimen is positive, the specimen is reported as “HBV DNA positive”.
- If an individual specimen is negative, the specimen is reported as “HBV DNA negative”.

Results of Pooled Source Plasma Specimens (Pools of up to 96 Individual Specimens)

The testing algorithm for testing of pooled specimens for the COBAS AmpliScreen HBV Test requires a single level of testing for Primary Pools that are negative for HBV DNA and three levels of testing (Primary Pool, Minipool, and confirmatory testing) for Primary Pools that are positive for HBV DNA.

Negative Primary Pools

When the Primary Pool is negative, report the results for all associated individual specimens in that Primary Pool as “HBV DNA Negative”.

Positive Primary Pools — Minipool Testing

Positive Primary Pools are traced to the positive individual using an overlapping pool testing matrix. Minipools are prepared from the eight individual specimens for columns 1 – 12 and from the 12 individual specimens for rows 1 – 8. The 20 minipools are tested using the Standard Specimen Processing Procedure. The positive specimen is identified by the intersection of the positive column and positive row. Confirmatory testing is conducted on the implicated specimen using the Standard Specimen Processing Procedure.

If the HBV DNA positive specimen is not identified because either the row or the column is negative, the original 96-specimen pool is tested in duplicate. The specimens are released as negative if the results of the duplicate testing are negative.

During deconstruction of a reactive pool,

1. If all subpools (whether independent or overlapping) are negative, all specimens are scored as non-reactive;
2. If deconstruction is done with non-overlapping subpools, the individual specimens of any non-reactive subpool are scored as non-reactive;

3. If deconstruction is done with overlapping pools, such as rows and columns, and if only one (or more) columns(s) is reactive, and no rows are positive, or vice versa, then all specimens in the reactive pool should be tested as individual specimens (i.e., this situation should not be resolved by retesting the master pool).

PROCEDURAL LIMITATIONS

1. This test has been evaluated only for use in combination with the COBAS AmpliScreen Multiprep Specimen Preparation and Control Kit, COBAS AMPLICOR Analyzer, and the Hamilton MICROLAB AT plus 2 Pipettor for the automated preparation of plasma pools.
2. ***Heparin inhibits PCR; specimens collected using heparin as the anticoagulant should not be used with the COBAS AmpliScreen HBV Test.***
3. Reliable results are dependent on adequate specimen collection and proper transport procedures.
4. Detection of HBV DNA is dependent on the number of virus particles present in the specimen and may be affected by specimen collection methods, patient factors (i.e., age, presence of symptoms), and/or stage of infection and pool size.
5. Only the Hamilton MICROLAB AT plus 2 Pipettor has been validated for use with the COBAS AmpliScreen HBV Test, for the automated preparation of plasma pools. Adhere to the hardware instructions and safety precautions outlined in the User Manual for the Hamilton MICROLAB AT plus 2 Pipettor.

PERFORMANCE CHARACTERISTICS

Reproducibility

The reproducibility of the COBAS AmpliScreen HBV Test was established by testing two 6-member EDTA plasma panels with known concentrations of HBV. Panel One, which was tested by using the Multiprep Specimen Processing Procedure, was composed

of HBV DNA positive specimens at concentrations of 25, 90, 150, and 25,000 copies/mL and two HBV negative specimens. Panel Two, which was tested by using the Standard Specimen Processing Procedure was composed of HBV positive specimens at concentrations of 75, 300, 500 and 25,000 copies/mL and two HBV negative specimens.

Testing was performed at three sites with two operators at each site using three COBAS AmpliScreen HBV Test kit lots and analyzed in 5 different days. Each operator used a dedicated COBAS AMPLICOR Analyzer throughout the study. Each operator was provided panel sets that had been randomized and labeled in blinded fashion.

All valid reproducibility data were evaluated by calculating the percentage of correct results for each panel member. The data were analyzed by site, lot, testing day, run, and operator for each Specimen Processing Procedure (Multiprep and Standard).

The reproducibility study for the COBAS AmpliScreen HBV Test demonstrated consistency by lot and site for both the Multiprep and Standard Specimen Processing Procedures as seen in Table 1 and Table 2 below. The reproducibility by operator is shown in Table 3 and Table 4 below.

Table 1: Reproducibility Results — Multiprep Specimen Processing Procedure

Results By Lot (# Positive / # Tested)					
	Negative	25 c/mL	90 c/mL	150 c/mL	25,000 c/mL
Lot #1	0/180	75/90	89/90	89/90	90/90
(%)	(0%)	(83%)	(99%)	(99%)	(100%)
Lot #2	0/178	75/90	87/88	89/90	90/90
(%)	(0%)	(83%)	(99%)	(99%)	(100%)
Lot #3	1/179	76/89	88/89	90/90	90/90
(%)	(1%)	(85%)	(99%)	(100%)	(100%)
Results By Site (# Positive / # Tested)					
Site #1	0/180	80/89	90/90	90/90	90/90
(%)	(0%)	(90%)	(100%)	(100%)	(100%)
Site #2	0/177	76/90	84/87	88/90	90/90
(%)	(0%)	(84%)	(97%)	(98%)	(100%)
Site #3	1/180	70/90	90/90	90/90	90/90
(%)	(1%)	(78%)	(100%)	(100%)	(100%)

Table 2: Reproducibility Results — Standard Specimen Processing Procedure

Results By Lot (# Positive / # Tested)					
	Negative	75 c/mL	300 c/mL	500 c/mL	25,000 c/mL
Lot #1 (%)	0/179 (0%)	76/89 (85%)	89/89 (100%)	90/90 (100%)	90/90 (100%)
Lot #2 (%)	1/179 (1%)	73/90 (81%)	88/89 (99%)	88/90 (98%)	90/90 (100%)
Lot #3 (%)	0/180 (0%)	78/90 (87%)	90/90 (100%)	90/90 (100%)	90/90 (100%)
Results By Site (# Positive / # Tested)					
Site #1 (%)	0/180 (0%)	72/89 (81%)	89/89 (100%)	90/90 (100%)	90/90 (100%)
Site #2 (%)	0/179 (1%)	76/90 (84%)	88/89 (99%)	88/90 (98%)	90/90 (100%)
Site #3 (%)	0/179 (0%)	79/90 (88%)	90/90 (100%)	90/90 (100%)	90/90 (100%)

Table 3: Operator Variability Data — Multiprep Specimen Processing Procedure

Results by Operator / Instrument (# Positive / # Tested)					
Operator	Negative	25 c/mL	90 c/mL	150 c/mL	25,000 c/mL
1	0/90 (0%)	39/44 (88.6%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
2	0/90 (0%)	41/45 (91.1%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
3	0/90 (0%)	44/45 (97.8%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
4	0/90 (0%)	32/45 (71.1%)	39/42 (92.2%)	43/45 (95.6%)	45/45 (100%)
5	0/90 (0%)	35/45 (77.8%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
6	1/90 (1.1%)	35/45 (77.8%)	45/45 (100%)	45/45 (100%)	45/45 (100%)

Table 4: Operator Variability Data — Standard Specimen Processing Procedure

Results by Operator / Instrument (# Positive / # Tested)					
Operator	Negative	75 c/mL	300 c/mL	500 c/mL	25,000 c/mL
1	0/90 (0%)	39/44 (88.6%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
2	0/90 (0%)	33/45 (73.3%)	44/44 (100%)	45/45 (100%)	45/45 (100%)
3	0/90 (0%)	41/45 (91.1%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
4	1/89 (1.1%)	35/45 (77.8%)	43/44 (97.7%)	43/45 (95.6%)	45/45 (100%)
5	0/89 (0%)	41/45 (91.1%)	45/45 (100%)	45/45 (100%)	45/45 (100%)
6	0/90 (0%)	38/45 (84.4%)	45/45 (100%)	45/45 (100%)	45/45 (100%)

Analytical Sensitivity — WHO HBV International Standard

The Limit of Detection of the COBAS AmpliScreen HBV Test was determined by using the WHO HBV International Standard (97/746). The WHO HBV International Standard was serially diluted in Negative Human Plasma to final concentrations of 100, 30, 10, 5, 4, and 3 IU/mL for the Multiprep Specimen Processing Procedure and 300, 100, 30, 20, 15, and 10 IU/mL for the Standard Specimen Processing Procedure. Each dilution was tested using two lots of COBAS AmpliScreen HBV Test at 60 replicates per lot.

When evaluated using PROBIT analysis, the combined data from all specimens tested with the Multiprep Specimen Processing Procedure indicate an average 95% LOD of 4.41 IU/mL, with lower and upper 95% confidence limits of 3.56 IU/mL and 6.13 IU/mL, respectively. The LOD of 4.41 IU/mL corresponds to approximately 22 copies/mL.

When evaluated using PROBIT analysis, the combined data from all specimens tested with the Standard Specimen Processing Procedure indicate an average 95% LOD of 15.99 IU/mL, with lower and upper 95% confidence limits of 13.78 IU/mL and 20.06 IU/mL, respectively. The LOD of 15.99 IU/mL corresponds to approximately 80 copies/mL.

Table 5 and Table 6 summarize the overall results for the Multiprep and Standard Specimen Processing Procedures, respectively.

Table 5: Serial Dilution Testing Summary for Multiprep Procedure with HBV DNA WHO International Standard (97/746) Combined Input Values with Lower 95% Confidence Limit (One-Sided)

HBV DNA Concentration (IU/mL)	Number of Positives	Number of Individual Tests	% Positives	95% Lower Confidence Limit One-Sided
100	120	120	100.0%	97.5%
30	120	120	100.0%	97.5%
10	120	120	100.0%	97.5%
5	115	120	95.8%	91.4%
4	108	120	90.0%	84.3%
3	112	120	93.3%	88.3%

Table 6: Serial Dilution Testing Summary for Standard Method with HBV DNA WHO International Standard (97/746) Combined Input Values with Lower 95% Confidence Limit (One-Sided)

HBV DNA Concentration (IU/mL)	Number of Positives	Number of Individual Tests	% Positives	95% Lower Confidence Limit One-Sided
300	120	120	100.0%	97.5%
100	120	120	100.0%	97.5%
30	118	119	99.2%	96.1%
20	116	120	96.7%	92.5%
15	115	120	95.8%	91.4%
10	101	120	84.2%	77.6%

Genotype Sensitivity and Inclusivity

One hundred fifteen specimens were diluted to 60 copies/mL and 200 copies/mL with HBV Negative Human Plasma (16 Genotype A, 22 Genotype B, 16 Genotype C, 8 Genotype D, 16 Genotype E, 22 Genotype F, 1 Genotype G, 3 Genotype A/C, 2 Genotype A/D, 2 Genotype C/D, 2 Genotype D/E, and 1 each of Genotypes A/E, B/C,

C/E, D/F, and E/F). All the genotypes tested positive by the COBAS AmpliScreen HBV Test at 60 copies/mL with the Multiprep Specimen Processing Procedure, and at 200 copies/mL with the Standard Specimen Processing Procedure. An additional 13 specimens (2 Genotype A, 2 Genotype B, 2 Genotype C, 2 Genotype D, 2 Genotype E and 3 Genotype F) were diluted in HBV Negative Human Plasma and tested in replicates of 22 with the COBAS AmpliScreen HBV Test at a level which resulted in a $\geq 95\%$ detection rate (2 to 100 copies/mL for the Multiprep Specimen Processing Procedure and 15 to 400 copies/mL for the Standard Specimen Processing Procedure). All of the genotypes tested positive by the COBAS AmpliScreen HBV Test. As a result of limited availability, only one Genotype G was evaluated. Therefore, performance of the COBAS AmpliScreen HBV Test with this one Genotype G specimen may not be representative of all Genotype G specimens. Data are provided in Table 7.

Table 7: HBV Genotypes Tested

Genotypes	Quantity	Reactive Total (Multiprep)	Reactive Total (Standard Prep)
A	18	18/18	18/18
B	24	24/24	24/24
C	18	18/18	18/18
D	10	10/10	10/10
E	18	18/18	18/18
F	25	25/25	25/25
G	1	1/1	1/1
A/C	3	3/3	3/3
A/D	2	2/2	2/2
C/D	2	2/2	2/2
D/E	2	2/2	2/2
A/E	1	1/1	1/1
B/C	1	1/1	1/1
C/E	1	1/1	1/1
D/F	1	1/1	1/1
E/F	1	1/1	1/1

Seroconversion Panels

Forty commercially available HBV seroconversion panels were tested with the COBAS AmpliScreen HBV Test. Each panel member was tested undiluted with the Standard Specimen Processing Procedure and diluted 1:24 with HBV Negative Human Plasma for testing with the Multiprep Specimen Processing Procedure. The COBAS AmpliScreen HBV Test results were compared to the results from an FDA-licensed HBV surface antigen (HBsAg) assay (Ortho HBsAg ELISA Test System 3).

When compared to the Ortho HBsAg ELISA Test System 3 assay, HBV DNA was consistently detected earlier than HBsAg by the Ortho Test in 38 of 40 panels using the Multiprep Specimen Processing Procedure (specimens diluted 1:24). With the Standard Specimen Processing Procedure (specimens undiluted), HBV DNA was consistently detected earlier than HBsAg by the Ortho Test in 39 of 40 panels. In no instance was there a specimen that showed HBsAg reactivity yet was negative for HBV DNA. The COBAS AmpliScreen HBV Test detected HBV DNA a mean of 15 ± 17 days prior to detection of HBsAg by the Ortho HBsAg test using the Multiprep Specimen Processing Procedure and 20 ± 17 days when using the Standard Specimen Processing procedure. See Table 8.

Table 8: Summary of the Pre-Seroconversion Detection of HBV DNA vs. HBsAg

	Days Prior to Ortho HBsAg Test System 3	
	Multiprep Specimen Processing Procedure (40 panels tested*)	Standard Specimen Processing Procedure (40 panels tested*)
Mean	15	20
Median	14	18
Maximum	30	53
Minimum	0	0

* One seroconversion panel was not included in the calculations which was detected by the COBAS AmpliScreen HBV Test 108+ days prior to the Ortho HBsAg ELISA Test System 3. In addition, the calculations do not include the results for five seroconversion panels tested with the Multiprep Specimen Processing Procedure and three seroconversion panels tested by the Standard Specimen Processing Procedure that were intermittently HBV DNA positive up to 100+ days prior to the Ortho test.

The results of these studies provide objective evidence that the COBAS AmpliScreen HBV Test, using either the Multiprep or Standard Specimen Processing Procedures, demonstrates greater sensitivity than observed with current HBsAg serology assays in detecting early HBV infection.

Analytical Specificity — Potentially Cross Reactive and Interfering Microorganisms

The analytical specificity of the COBAS AmpliScreen HBV Test was evaluated by testing a panel of 38 microorganisms, including 33 viral isolates and 5 bacterial strains, using the Multiprep Specimen Processing Procedure. No cross-reactivity was observed with the COBAS AmpliScreen HBV Test. Table 9 provides a list of organisms tested.

Table 9: Analytical Specificity — Micororganisms Tested

Adenovirus, Human Type 2	Chlamydia trachomatis	HCV 1a
Adenovirus, Human Type 3	Neisseria gonorrhoeae	HCV 1b
Adenovirus, Human Type 7	Epstein Barr Virus (Burkitt's lymphoma)	HCV 2a/2c
Cytomegalovirus (3 strains)	Epstein Barr Virus (RAJI Human Burkitt's lymphoma)	HCV 2b
Herpes Simplex type 1	Echovirus 1	HCV 3a
Herpes Simplex type 2	Echovirus 5	HIV Subtype A
Hepatitis A	Coxsackievirus B1	HIV Subtype B
Human Papilloma Virus, Type 6a	Varicella-Zoster	HIV Subtype C
Human Papilloma Virus, Type 16	Varicella	HIV Subtype D
Human Papilloma Virus, Type 18	Propionibacterium acnes	HIV Subtype E
HTLV-I	Staphylococcus aureus	HIV Subtype F
HTLV-II	Staphylococcus epidermidis	HIV Subtype G

Up to 25 individual patient plasma specimens from each of the following disease categories were spiked with HBV positive plasma to a level of 90 copies/mL for the Multiprep Specimen Processing Procedure and 300 copies/mL for the Standard Specimen Processing Procedure or non-spiked: specimens that are serologically positive for HIV-1, HIV-2, HAV, HCV, autoimmune disease, EBV, CMV, and negative human plasma spiked with Candida albicans. No false negative test results were observed with the HBV

spiked specimens, and no cross-reactivity was observed with the non-spiked specimens, using both the Multiprep and Standard Specimen Processing Procedures.

Potentially Interfering Substances

Endogenous Interfering Substances

HBV spiked and non-spiked plasma samples derived from whole blood containing abnormally high concentrations of bilirubin (up to 20 mg/mL), triglycerides (up to 3000 mg/dL), hemoglobin (up to 5.0 g/dL), and albumin (up to 6 g/dL) were tested. These endogenous substances did not interfere with the sensitivity or specificity of the COBAS AmpliScreen HBV Test using either the Multiprep or Standard Specimen Processing Procedure.

Exogenous Interfering Substances

HBV spiked and non-spiked plasma samples derived from whole blood containing abnormally high concentrations of aspirin (up to 50 mg/mL), pseudoephedrine-HCl (up to 3 mg/dL), ascorbic acid (up to 20 mg/dL), acetaminophen (up to 40 mg/dL), or ibuprofen (up to 40 mg/dL) were tested. These exogenous substances did not interfere with the sensitivity or specificity of the COBAS AmpliScreen HBV Test using either the Multiprep or Standard Specimen Processing Procedures.

CLINICAL PERFORMANCE

Detection of HBV DNA in Seropositive Specimens

A total of 1177 known HBV seropositive (HBsAg positive) specimens were tested by the COBAS AmpliScreen HBV Test. These HBV seropositive specimens included the following specimens and sources: 723 HBV seropositive specimens, including 49 acute patients (HBsAg and HBeAg positive) obtained from commercial vendors and blood banks in the US, Japan and China; 100 chronic HBV patient specimens (HBsAg positive for at least 6 months) obtained from a commercial vendor; 60 HBsAg seropositive

specimens collected from patients at high risk for hepatitis, and 189 HBsAg positive specimens from 40 seroconversion panels (also described in Seroconversion Panels). These specimens were tested using both the Multiprep Specimen Processing Procedure (specimens diluted 1:24 in negative human plasma) and Standard Specimen Processing Procedure (specimens tested undiluted). An additional 105 HBsAg positive specimens were tested in the Clinical Performance study. These specimens were initially tested in primary mini-pools of 24 specimens using the Multiprep Specimen Processing Procedure, and tested individually using the Standard Specimen Processing Procedure.

Seropositive Specimens (Including 49 Acute Patients) from Commercial Vendors

A total of 723 HBV seropositive specimens, including 49 acute patients (HBsAg and HBeAg positive) obtained from commercial vendors and blood banks in the US, Japan, and China were tested with the COBAS AmpliScreen HBV Test.

Of the 723 HBsAg positive specimens that were tested using the Multiprep Specimen Processing Procedure, 694 specimens (96.0%) tested HBV DNA positive and 29 specimens (4.0%) were HBV DNA negative. Of the 723 HBsAg positive specimens that were tested individually using the Standard Specimen Processing Procedure, 708 specimens (97.9%) tested HBV DNA positive and 15 specimens (2.1%) tested HBV DNA negative. The specimens discordant with serology when tested individually using the Standard Specimen Processing Procedure were a subset of the specimens which were discordant with serology when tested using the Multiprep Specimen Processing Procedure. *See* Table 10 and Table 11.

Table 10: HBV Seropositive Specimens (Including 49 Acute Patients) from Commercial Vendors —Results for Multiprep Specimen Processing Procedure (Specimens Diluted 1:24)

		HBV DNA Results
Multiprep Procedure	+	694 (96.0%)
	–	29*
Total		723

* 22 of the 29 HBV DNA negative specimens were tested by alternate HBV DNA tests and were negative (< 300 copies/mL). Includes 7 Specimens which were not retested by alternate HBV DNA tests due to insufficient volume.

Table 11: HBV Seropositive Specimens Including 49 Acute Patients from Commercial Vendors —Results for Standard Specimen Processing Procedure (Specimens Undiluted)

		HBV DNA Results
Standard Procedure	+	708 (97.9%)
	–	15*
Total		723

* All HBV DNA negative specimens were tested by alternate HBV DNA tests and contained < 300 copies/mL HBV DNA.

Chronic HBV Patients

A total of 100 HBV seropositive specimens from chronic HBV patients were tested with the COBAS AmpliScreen HBV Test. Of the 100 HBsAg positive specimens, there were 78 specimens (78.0%) HBV DNA positive and 22 specimens (22.0%) HBV DNA negative with the Multiprep Specimen Processing Procedure, and 95 specimens (95.0%) HBV DNA positive and 5 specimens (5.0%) HBV DNA negative with the Standard Specimen Processing Procedure. The specimens discordant with serology when tested individually using the Standard Specimen Processing Procedure were a subset of the specimens which were discordant with serology when tested using the Multiprep Specimen Processing Procedure. See Table 12 and Table 13.

**Table 12: HBV Seropositive Specimens from Chronic HBV Patients —
Results for Multiprep Specimen Processing Procedure
(Specimens Diluted 1:24)**

		HBV DNA Results
Multiprep Procedure	+	78 (78.0%)
	–	22*
Total		100

* All HBV DNA negative specimens were tested by alternate HBV DNA tests and contained < 300 copies/mL HBV DNA

**Table 13: HBV Seropositive Specimens from Chronic HBV Patients —
Results for Standard Specimen Processing Procedure
(Specimens Undiluted)**

		HBV DNA Results
Standard Procedure	+	95 (97.9%)
	–	5*
Total		100

* All HBV DNA negative specimens were tested by alternate HBV DNA tests and and contained < 300 copies/mL HBV DNA.

High Risk Patients

There were 60 HBV seropositive specimens from patients at high risk for liver disease tested with the COBAS AmpliScreen HBV Test. There were 59 specimens (98.3%) which tested HBV DNA positive and 1 specimen (1.7%) HBV DNA negative using the Multiprep Specimen Processing Procedure. All 60 specimens (100%) were HBV DNA positive using the Standard Specimen Processing Procedure. *See Table 14 and Table 15.*

Table 14: HBV Seropositive Specimens from Patients at High Risk for Liver Disease —Results for Multiprep Specimen Processing Procedure (Specimens Diluted 1:24)

		HBV DNA Results
Multiprep Procedure	+	59 (98.3%)
	–	1*
Total		60

* The HBV DNA negative specimen was tested by alternate HBV DNA tests and contained < 300 copies/m HBV DNA

Table 15: HBV Seropositive Specimens from Patients at High Risk for Liver Disease —Results for Standard Specimen Processing Procedure (Specimens Undiluted)

		HBV DNA Results
Standard Procedure	+	60 (100%)
	–	0
Total		60

Seropositive Specimens from Seroconversion Panels

For the 40 HBV seroconversion panels, each panel was consistently positive by the COBAS AmpliScreen HBV Test with both the Multiprep and Standard Specimen Processing Procedures. Combining all 40 panels, a total of 189 HBsAg positive specimens were tested and shown to be 100% concordant (HBV DNA positive) with serologic results.

Seropositive Specimens from Blood Donors

During the clinical studies, there were 105 HBV specimens confirmed positive for HBsAg which were also tested with the COBAS AmpliScreen HBV Test. There were 87 specimens (82.9%) HBV DNA positive and 18 specimens (17.1%) HBV DNA negative using the Multiprep Specimen Processing Procedure. Of these 105 HBsAg

positive specimens, 104 specimens were tested individually using the Standard Specimen Processing Procedure, and 97 specimens (93.3%) tested DNA positive and 7 specimens (6.7%) tested DNA negative. The specimens discordant with serology when tested individually using the Standard Specimen Processing Procedure were a subset of the specimens which were discordant with serology when tested using the Multiprep Specimen Processing Procedure. *See* Table 16 and Table 17.

Table 16: HBV Seropositive Specimens — Results for Multiprep Specimen Processing Procedure (Specimens Diluted 1:24)

		HBV DNA Results
Multiprep Procedure	+	87 (82.9%)
	–	18*
Total		105

** 6 of 18 specimens were negative for HBV DNA by alternate qualitative NAT and specimens from 3 of 18 were not available for testing by alternate qualitative NAT. 11 of 18 were quantitated for HBV DNA by alternate quantitative NAT and 7 of these contained < 300 copies/mL HBV DNA.*

Table 17: HBV Seropositive Specimens — Results for Standard Specimen Processing Procedure (Specimens Undiluted)

		HBV DNA Results
Standard Procedure	+	97 (93.3%)
	–	7*
Total		104

** 4 of 7 specimens were negative for HBV DNA by alternate NAT 1 was not tested by alternate NAT and the remaining 2 quantified by alternate NAT had HBV DNA < 100 copies/mL).*

Pool Reactivity in Whole Blood

A random selection of 25,845 pools revealed that 150 Primary Pools were reactive for an initial reactive rate of 0.58%. There were 85/150 (56.7%) positive pools that were concordant with confirmed positive serology status. Two of these pools were identified as having a window period case. There were 51/25,845 (0.197%) pools that were initially reactive but determined to be HBV DNA negative upon resolution testing. A total of

9/25,845 (0.035%) pools were found positive but were not confirmed positive by serology or by subsequent testing of individual specimens by the COBAS AmpliScreen HBV Test. There were 3 pools with two positive units, one with concordant serology and one without concordant serology. Results are summarized in Table 18.

Table 18: Pool Reactivity in Whole Blood

Category	No. of Pools	Percentage
Pools tested	25,845	100%
Non-Reactive pools	25,695	99.42%
Initially Reactive pools	150	0.58%
Initially reactive pools with positive individual specimen testing and with concordant serology	85	0.33%
Positive pools due to window period case	2	0.008%
Initially reactive pools with negative serology and negative individual specimen AmpliScreen Testing (false positive)	51	0.2%
Initially reactive pools with positive individual specimen testing and without concordant serology	9	0.03%
Initially reactive pools with 2 positive units by individual testing; one concordant with serology and one without concordant serology	3	0.01%

Of the 25,845 pools, a total of 581,790 specimens were tested from 5 geographically divergent sites. The results from these specimens were used to determine the specificity and sensitivity of COBAS AmpliScreen HBV Test. The HBV serology status of each specimen was determined using each specimens antigen and antibody results. HBV serology status-negative included specimens that were HBsAg and anti-HBc negative unless the subject was enrolled in the follow-up study and had test results that changed this assessment. HBV serology status positive included those specimens that were HBsAg positive regardless of the anti-HBc results unless the subject was enrolled in the follow-up study and had test results that changed this assessment. HBV serology status unknown included those specimens that were anti-HBc positive and HBsAg negative.

There were 578,694 specimens that were determined to be HBV serology status negative (2 HBV MP NAT positive serology negative). Of these, 578,673 were also HBV DNA negative (21 specimens were false positive). The specificity of the COBAS AmpliScreen HBV Test in this study was 578,673/578,694 or 99.9964% with 95% confidence limits of 99.9948% to 99.9979%.

There were 105 specimens that were determined to be HBV serology status positive (HBsAg confirmed positive). Of these 105 HBsAg specimens, 87 specimens (82.9%) were also HBV DNA positive when tested in primary mini-pools using the Multiprep Specimen Processing Procedure and 18 specimens (17.1%) were negative by this procedure. Of these 18 HBV DNA negative specimens, 17 specimens were also tested individually by the Standard Specimen Processing Procedure and 10 of these specimens were positive by this procedure and 7 specimens were negative.

It should be noted that three of the specimens with quantitative detection values of 1200 copies/ml, 2600 copies/ml and 5900 copies/ml of DNA that were also positive for HBsAg and anti-HBc were negative on mini-pool NAT.

Detection of Window Period Cases

A confirmed window period case is defined as an individual from whom the index donation was positive in the COBAS AmpliScreen HBV Test but tested negative by HBsAg and anti-HBc and a follow-up specimen tested positive either by the COBAS AmpliScreen HBV Test, HBsAg, or anti-HBc. Two window period cases were found during the clinical trial for a detection rate of 1:290,895 rate with exact 95% confidence limits of 1:11,497,826 to 1:52,083.

Single Donation Testing Performance

A total of 1754 specimens for which serology results were available were tested individually in the COBAS AmpliScreen HBV Test clinical trial. The results are shown in Table 19. There were 89/1754 classified as HBV status positive (87 were HBsAg, anti-HBc and HBV DNA positive and 2 were HBV window period cases based on follow-up testing). In the absence of follow-up testing 13/1754 were classified as HBV status unknown (12 were HBsAg negative, anti-HBc positive and HBV DNA negative, and 1 was HBsAg negative, anti-HBc positive and HBV DNA positive). There were 1652/1754 classified as HBV status negative and of these, 1625/1652 were COBAS AmpliScreen HBV Test negative. Of the remaining 27 specimens that tested positive by the COBAS AmpliScreen HBV Test and negative by serology testing, 12 were enrolled in the follow-up study and determined to be false positives. Based on these results, all 27 were presumed to be false positives on the COBAS AmpliScreen HBV Test. The specificity of the COBAS AmpliScreen HBV Test in this study was 98.4% (1625/1652) with a 95% confidence interval of 97.6% to 98.9%.

Table 19: Paired Specimen Results for Individual Samples with Assigned HBV Status

HBsAg Result	Anti-HBc Result	COBAS AmpliScreen DNA Result	Status	Total
Negative	NR	Negative	Negative	1625
Negative	RR	Negative	Unknown	12
Negative	NR	Positive	Negative	27
Negative	NR	Positive	Positive*	2
Negative	RR	Positive	Unknown	1
Positive	Any	Negative	Positive	0
Positive	Any	Positive	Positive	87
			Total	1754

* Status Positive reclassified due to follow-up results.

Performance Characteristics for Source Plasma

Clinical Performance

A total of 1,080 primary pools tested in the 96-member mini-pool format representing 103,680 specimens from 40,230 donors revealed that 8 pools were reactive with the COBAS AmpliScreen HBV Test for an initial reactive rate of 0.74%. Of the 8 reactive pools, there were 3 identified HBV DNA positive pools and 2 pools were positive due to apparent window period cases. The remaining 3 pools were reactive for HBV DNA, but were not confirmed. The data are presented in Table 20.

Table 20: Pool Reactivity in Source Plasma

Category	No. of Pools	Percentage
Pools tested	1080	100
Non-reactive pools	1072	99.26
Initially reactive pools	8	0.74
Initial pools containing a reactive individual specimen with concordant serology	3	0.28
Positive pools due to window period case ¹	2	0.18%
Initially reactive pools with negative resolution COBAS AmpliScreen Testing (false positive)	3	0.28%

¹ Two HBsAg negative specimens were in one 96-member mini-pool

There were 1075 pools used to determine the specificity of the COBAS AmpliScreen HBV Test. Of these pools, 1072 were HBV DNA negative and 3 were Initially Reactive with negative resolution COBAS AmpliScreen Testing (false positive). The specificity of the COBAS AmpliScreen HBV Test in this study was 1072/1075 or 99.7209% with 95% confidence limits of 99.19% to 99.94%.

HBV Seroconversion Panels

Ten commercially available HBV seroconversion panels were tested using the Multiprep Specimen Processing Procedure and compared to results obtained with licensed HBsAg tests. Blinded panel members were diluted 1:96 with HBV negative human plasma.

In two panels, COBAS AmpliScreen HBV Test detected HBV DNA on the same day as HBsAg was detected by the Ortho Antibody to HBsAg ELISA Test System 2. In three panels, COBAS AmpliScreen HBV Test detected HBV DNA on the same day as HBsAg was detected by the Ortho HBsAg Test System 3. In two panels, COBAS AmpliScreen HBV Test detected HBV DNA on the same day as HBsAg was detected by the Abbott Auszyme Test. Data are presented in Table 21.

Table 21: Summary of Pre-Seroconversion Detection of HBV DNA vs. FDA Licensed Serology Tests — Multiprep Specimen Processing Procedure (Specimen diluted 1:96)

	Days Before Ortho HBsAg ELISA Test System 2 (8 panels tested*)	Days BeforeOrtho HBsAg ELISA Test System 3 (9 panels tested*)	Days Before Abbott Auszyme overnight for HBsAg (9 panels tested*)
Mean	8.8	8.3	11.0
Median	8	7	9
Maximum	27	27	27
Minimum	0	0	0

* One seroconversion panel was not included in the calculations which was detected by the COBAS AmpliScreen HBV Test 108+ days prior to detection of HBsAg by the Ortho HBsAg ELISA Test System 2, 101 days prior to the Ortho HBsAg ELISA Test System 3, and 87 days prior to the Abbott Auszyme test.

The seroconversion study demonstrates the COBAS AmpliScreen HBV test used with the Multiprep Specimen Processing Procedure and pools of 96 specimens, identified HBV infected specimens at the same time or earlier than the U.S. FDA licensed HBsAg tests.

REFERENCES

1. Fagan, EA and Harrison, TJ. 2000. Viral Hepatitis: A Handbook for Clinicians and Scientists, Springer-Verlag New York Inc., 89-130.
2. Beasley RP. 1982. Hepatitis B virus as etiologic agent in hepatocellular carcinoma – epidemiological considerations. *Hepatology* 2 (suppl):21S-26S.
3. Feitelson M. 1992. Hepatitis B virus infection and primary hepatocellular carcinoma. *Clin Microbiol Rev* 14:257-301
4. Mast EE and Later MJ. 1993. Epidemiology of viral hepatitis: an overview. *Semin Virol* 4:273-283
5. Robinson WS, Clayton DA, and Greenman RL. 1974. DNA of a human hepatitis B virus candidate. *J Virol* 14:384-391
6. Galibert F, Mandart E, Fitoussi F, Tiollais P, and Charney P. 1979. Nucleotide sequence of the hepatitis B virus genome (subtype ayw) cloned in *E. coli*. *Nature* 281:646-650
7. Kitchen A. 1998. Hepatitis B and blood safety. *Vaccine* 16:s34-s37.
8. Schreiber GB, Busch MP, Kleinman SH, and Korelitz JJ. 1996. The risk of transfusion-transmitted viral infections. *N Engl J Med* 334:1685-1690
9. Zuckerman AJ. 1999. More than third of the world's population has been infected with hepatitis B virus [letter]. *Br Med J* 318:1213
10. McQuillan GM, Coleman PJ, Kruszon-Moran D, Moyer LA, Lambert SB, and Margolis HS. 1999. Prevalence of hepatitis B virus infection in the United States: The National Health and Nutrition Examination Survey, 1976 through 1994. *Am J Pub Health* 89:14-18
11. Van DP and Vellinga A. Epidemiology of hepatitis B and C in Europe. *Acta Gastroenterol. Belg.* 1998;61:175-182

-
12. Maddrey WC. 2000. Hepatitis B: An important public health issue. *J Med Virol* 61:362-366
 13. Larsen J, Hetland G, and Skaug K. 1990. Posttransfusion hepatitis B transmitted by blood from a hepatitis B surface antigen-negative hepatitis B virus carrier. *Transfusion* 30:431
 14. Saraswat S, Banerjee K, Chaodhury N, Mahant T, Khandekar P, Gupta RK, and Naik S. 1996. Post-transfusion hepatitis type B following multiple transfusions of HBsAg-negative blood. *J of Hepatol* 25:639-643
 15. Yotsuyanagi H, Yasuda K, Iino S, Moriya K, Shintani Y, Fujie H, Tsutsumi T, Kimura S, and Koike K. 1998. Persistent Viremia After Recovery From Self-Limited Acute Hepatitis B. *Hepatology* 27:1377-1382
 16. Michalak TI, Pasquinelli C, Guilhot S, and Chisari FV. 1994. Hepatitis B virus persistence after recovery from acute viral hepatitis. *J Clin Invest* 93:230-239
 17. Cabrerizo M, Bartolomé J, Caramelo C, Barril G, AND Carreño V. 2000. Molecular Analysis of Hepatitis B Virus DNA in Serum and Peripheral Blood Mononuclear Cells From Hepatitis B Surface Antigen–Negative Cases. *Hepatology* 32:116-123
 18. Busch MP, Stramer SL, Kleinman SH. 1997. Evolving applications of nucleic acid amplification assays for prevention of virus transmission by blood components and derivatives. In: Garratty G, ed. *Applications of molecular biology*. Bethesda, American Association of Blood Banks 1997:121-73. (presented at a workshop during the 50th Annual Meeting of the AABB, October 1997, Denver, CO)
 19. Mortimer J. 1997. Intersecting pools and their potential application in testing donated blood for viral genomes. *Vox Sanguinis* 73:93-96

-
20. Kanemitsu K, Tomono T, Murozuka T, Emura H, Yugi H, Miyamoto, and Nishioka K. 2000. National NAT screening of HBV, HCV, and HIV for blood transfusion. *Vox Sanguinis* 78 S1:0081
 21. Yang Y, Xu D, Mendoza M, Lamendola M, Wu Y, Yeh S, Kung K. 2000. A prototype of Hepatitis B virus (HBV) PCR assay with a semi-automated instrument system. *Transfusion* 40 No.10S:SP207
 22. Rosenstraus M, Gutekunst K, Herman S et al. 1998. Improved COBAS AMPLICOR® viral assays as a basis for mini-pool screening of viruses in blood or plasma. *Infusionsther Tranfusionsmed* 25:153-9
 23. Longo M.C., Berninger, M.S. and Hartley, J.L. 1990. Use of uracil DNA glycosylase to control carry-over contamination in polymerase chain reactions. *Gene* 93:125-128.
 24. Richmond, J.Y. and McKinney, R.W. eds. 1999. Biosafety in Microbiological and Biomedical Laboratories. HHS Publication Number (CDC) 93-8395.
 25. National Committee for Clinical Laboratory Standards. Protection of Laboratory Workers from Infectious Disease Transmitted by Blood, Body Fluids, and Tissue. Approved Guideline. NCCLS Document M29-A Villanova, PA:NCCLS, 1997.
 26. International Air Transport Association. Dangerous Goods Regulations, 41st Edition. 2000. 704 pp.
 27. Centers for Disease Control and Prevention. Summary of Notifiable Diseases, United States, 2000, Morbidity and Mortality Weekly Report, 2002;49(53).

Manufactured by:

Roche Molecular Systems, Inc., Branchburg, NJ 08876 USA

A Member of the Roche Group

Distributed by:

Roche Diagnostics Corporation
Indianapolis, IN 46256 USA
(For Technical Assistance call the
Roche Response Center®
toll-free: 1-800 526 1247)

Roche Diagnostics
201 Boul. Armand-Frappier
Laval, Quebec, Canada H7V 4A2
(For Technical Assistance call:
Pour toute assistance technique,
appeler le: 1-800 268 0482)

Roche Diagnostics (Schweiz) AG
CH-6343 Rotkreuz
Roche Diagnostics GmbH
D-68298 Mannheim, Germany

Roche Diagnostics
F-38240 Meylan
Distributore in Italia
Roche Diagnostics SpA
Piazza Durante 11
I-20131 Milano

Roche Diagnostics S.L.
E-08006 Barcelona

Distribuidor em Portugal:
Roche Farmacêutica Química, Lda
P-2700 Amadora

United States:

Roche® is a registered trademark of Hoffmann-La Roche Inc.

AmpErase® is a registered trademark and COBAS Amplicor™, COBAS AmpliScreen™ and AMPLILINK™ are trademarks of Roche Molecular Systems, Inc., licensed to Roche Diagnostics Corporation.

Roche Response Center® is a registered service mark of Roche Diagnostics Corporation.

Other Countries:

The above trademarks are either owned by or used under license by the Roche affiliate distributing this product.

Dynabeads® paramagnetic particles are licensed under patents owned by Dynal Biotech ASA, Oslo, Norway. Dynabeads® is a registered trademark of Dynal Biotech ASA, Oslo, Norway, licensed to Roche Diagnostics Corporation, Indianapolis, Indiana.

Eppendorf®, Eppendorf Multipette® and Eppendorf Combitip® are registered trademarks of Eppendorf-Netheler-Hinz GmbH, Hamburg, Germany.

MicroAmp® is a registered trademark of Applied Biosystems Corporation or its subsidiaries in the U.S. and certain other countries

ProClin® is a registered trademark of Rohm and Haas Company.

Copyright 200X, Roche Molecular Systems, Inc. All rights reserved.